
CE
0120

COVERED CP STENT™

COARCTATION OF THE AORTA

Instructions for Use.....	2
Mode d'emploi	7
Istruzioni per l'uso.....	12
Gebrauchsanweisung.....	18
Instrucciones de utilización.....	24
Bruksanvisning	29
Brugsanvisning	34
Gebruiksaanwijzing	39
Instruções de utilização	45
Kullanım Talimatı	51

CAUTION: Federal (USA) Law restricts this device to sale by or on the order of a physician.

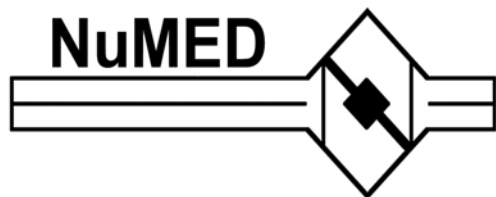Instructions for Use:

INDICATIONS

Indicated for implantation in the native and/or recurrent coarctation of the aorta on patients with the following clinical conditions:

- Stenosis of the aorta resulting in significant anatomic narrowing as determined by angiography or non-invasive imaging, i.e. echocardiography, magnetic resonance imaging (MRI), CT scan;
- Stenosis of the aorta resulting in hemodynamic alterations, resulting in systolic pressure gradient, systemic hypertension or altered left ventricular function;
- Stenosis of the aorta where balloon angioplasty is ineffective or contraindicated;
- Stenosis diameter <20% of adjacent vessel diameter. Stenosis that would present increased risk of vascular damage or disruption; or aneurysm associated with coarctation of the aorta.

DESCRIPTION

The Covered CP Stent is balloon expandable and intended for permanent implant. The Covered CP Stent is composed of heat treated 90% platinum/10% iridium wire that is arranged in laser welded rows with a "zig" pattern. The number of zigs in a row can be varied and will impact the strength of the stent as well as the eventual expanded diameter and percent stent shortening, while the number of rows will determine the unexpanded length of the stent. The Covered CP Stent has an ePTFE covering attached to the stent framework. This covering acts as a fluid barrier creating a fluid tight conduit through the stent length.

HOW SUPPLIED

Supplied sterilized by ethylene oxide gas. Sterile and non-pyrogenic if package is unopened or undamaged. Do not use the product if there is doubt as to whether the product is sterile. Avoid extended exposure to light. Upon removal from package, inspect the product to ensure no damage has occurred.

CONTRAINDICATIONS

- Patients too small to allow safe delivery of the stent without compromise to the systemic artery used for delivery;
- Unfavorable aortic anatomy that does not dilate with high pressure balloon angioplasty;
- Occlusion or obstruction of systemic artery precluding delivery of the stent;
- Clinical or biological signs of infection;
- Active endocarditis;
- Known allergy to aspirin, other antiplatelet agents, or heparin;
- Pregnancy.

WARNINGS

- As with any type of implant, infection secondary to contamination of the stent may lead to aortitis, or abscess. The platinum/iridium stent may migrate from the site of implant. Over-stretching of the artery may result in rupture or aneurysm formation.
- When the stent is crimped onto a balloon delivery catheter, the maximum balloon inflation pressure must not exceed the recommended inflation pressure specified in the manufacturer's instructions.
- The inflation diameter of the balloon used during stent delivery should approximate the diameter of the obstructive vessel and the intended implant site.
- Excessive force while crimping may weaken welds of the stent.
- This device is intended for single use only. Do not resterilize and/or reuse it, as this can potentially result in compromised device performance and increased risk of cross-contamination.

PRECAUTIONS

- Use of an inflation device with pressure gauge is highly recommended during this procedure.
- The stent is rigid and may make negotiation through vessels difficult.
- Dilatation procedures should be conducted under fluoroscopic guidance with appropriate x-ray equipment.
- Guidewires are delicate instruments. Care should be exercised while handling to help prevent the possibility of breakage.
- Careful attention must be paid to the maintenance of tight catheter connections and by aspiration before proceeding to avoid air introduction into the system.
- Under no circumstances should any portion of the catheter system be advanced against resistance. The cause of the resistance should be identified with fluoroscopy and action taken to remedy the problem.

POTENTIAL COMPLICATIONS/ADVERSE EFFECTS

NOTE: Circumferential tear of the delivery balloon catheter prior to complete expansion of the stent may cause the balloon to become tethered to the stent, requiring surgical removal. In case of rupture of an adequately sized balloon after stent expansion, it can be withdrawn and a new balloon catheter exchanged over a guidewire to complete expansion of the stent.

Cardiac catheterization carries certain risks. In addition, potential complications, and related adverse effects associated with implants include, but are not limited to:

- Stent migration
- Stent Fracture
- Vessel ruptures
- Minor hematoma
- Distal thromboemboli
- Death
- Intraluminal thrombus
- Hematoma requiring repair
- Bleeding requiring transfusion
- Cell necrosis at the site of implant
- Stent misplacement
- Pseudoaneurysm
- Stent misorientation
- Sepsis/infection
- AV fistula formation
- Transitory arrhythmia
- Cerebrovascular incident

MRI SAFETY INFORMATION

Nonclinical testing and modeling has demonstrated that the CP Stent is MR Conditional. A patient with this device can be safely scanned in an MR system meeting the following conditions:

- Static magnetic field of 1.5 T and 3 T
- Maximum spatial gradient magnetic field of 2500 gauss/cm (25 T/m)
- Maximum MR system reported, whole body averaged specific absorption rate (SAR) of 2.0 W/kg for 15 minutes of scanning (Normal Operating Mode)

Based on nonclinical testing and modeling, under the scan conditions defined above, the CP Stent is expected to produce a maximum in vivo temperature rise of less than 2°C after 15 minutes of continuous scanning.

MR image quality may be compromised if the area of interest is in the same area, or relatively close to the position of the device. In nonclinical testing, the image artifact caused by the device extends approximately 3 mm from the CP Stent when imaged with a spin echo pulse sequence and 6 mm when imaged with a gradient echo pulse sequence and a 3 T MRI System. The lumen of the device was obscured.

The presence of other implants or medical circumstances of the patient may require lower limits on some or all of the above parameters.

INSTRUCTIONS FOR USE

Select Stent Size

1. Measure the length of the target stricture to determine the length of stent required. Size the stent length to extend slightly proximal and distal to the stricture.
2. The appropriate stent length should be selected based on covering the entire obstructed segment with a single stent.
Note: Should more than one stent be required, place the stent most distal from the puncture site first, followed by placement of the proximal stent in tandem.
3. Measure the diameter of the reference stricture and vessel proximal and distal to the target lesion to determine the appropriate size stent and delivery system.

Preparation of Stent

- Remove the stent and rinse in sterile saline.

Preparation of Stent Delivery Catheter

- See instruction for use provided with the recommended balloon catheter.

Preparation of Stent Delivery System

1. Always place an appropriate sized guidewire through the balloon catheter lumen while prepping the balloon and crimping the stent.
2. Remove the balloon protector and visually inspect the balloon being used to insure that it is properly folded to its lowest profile in preparation for the stent. "Dry Prepping" the balloon delivery catheter by using negative pressure is highly recommended.
3. Slide the stent over the distal end of the balloon, maintaining the balloon fold, until the radiopaque markers are equal distance from the ends of the stent. Check the stent position under the fluoroscope before crimping.
4. Gently crimp (see insert IFU-CPCE) the stent on the balloon by using finger pressure and a "rolling action" to exert equal pressure on all sides of the stent. Only crimp until no movement is felt on the catheter. Platinum is very malleable and should be easily crimped. Place a small amount of undiluted contrast to "coat" the stent and improve adherence to the balloon. (AVOID BENDING OR TWISTING THE STENT).
5. Visually inspect the balloon/stent assembly to assure proper placement of the stent (MAY BE CONFIRMED BY FLUOROSCOPY).

Stent Deployment

1. The assembly is advanced through the long delivery sheath and over the stiff guidewire into the desired location for implant.
2. After correct positioning of the stent, pull back on the sheath to expose the stent. Confirm proper stent position by a small injection of contrast through the sidearm of sheath or through a second catheter. Expand the stent by inflating the balloon catheter to the recommended inflation pressure. Do not exceed the balloon manufacturer's rated burst pressure.

Delivery System Withdrawal

1. Once the stent is expanded, deflate the balloon and rotate to insure the stent is free and properly deployed.
2. Remove the balloon catheter and confirm the result with angiography.

NOTE: Diameter of the stent may be increased after placement by expanding with a larger diameter balloon. Do not exceed the maximum recommended expanded stent diameter.

RETURN OF EXPLANTED DEVICE:

NuMED, Inc. is interested in obtaining recovered CP Stents. Place the explanted device in a container or vial immediately after excision. For further instructions on the return of an explanted device, contact RA Manager, NuMED, Inc. 2880 Main Street, Hopkinton, New York, 12965. Phone number: 315-328-4491.

WARNING:

NuMED stents are placed in the extremely hostile environment of the human body. Stents may fail to function for a variety of causes including, but not limited to, medical complications or failure of stent by fracture and embolization. In addition, despite the exercise of all due

care in design, component selection, manufacture, and testing prior to sale, stents may be easily damaged before, during, or after insertion by improper handling, crimping or other intervening acts. Metal stents placed where there are extrinsic forces of compression, i.e. right ventricular outflow tract, are especially prone to fatigue fracture and embolization and should be avoided.

Warranty and Limitations

Stents and accessories are sold in an ‘as is’ condition. The entire risk as to the quality and performance of the stent is with the buyer. NuMED disclaims all warranties, expressed or implied, with respect to catheters and accessories, including but not limited to, any implied warranty of merchantability or fitness for a particular purpose. NuMED shall not be liable to any person for any medical expenses or any direct or consequential damages resulting from the use of any catheter or accessory or caused by any defect, failure, or malfunction of any catheter or accessory, whether a claim for such damages is based upon warranty, contract, tort, or otherwise. No person has any authority to bind NuMED to any representation or warranty with respect to catheters and accessories.

CP Stent™ Foreshortening Chart

Inflated Balloon Diameter	CP8Z16 (Stent length after expansion) Percentage Shortening	CP8Z22 (Stent length after expansion) Percentage Shortening	CP8Z28 (Stent length after expansion) Percentage Shortening	CP8Z34 (Stent length after expansion) Percentage Shortening	CP8Z39 (Stent length after expansion) Percentage Shortening	CP8Z45 (Stent length after expansion) Percentage Shortening
12mm	(1.61) cm 2.8%	(2.18) cm 0.8%	(2.62) cm 4.4%	(3.23) cm 3.1%	(3.72) cm 1.9%	(4.17) cm 3.8%
14mm	(1.54) cm 6.5%	(2.08) cm 5.4%	(2.56) cm 6.8%	(3.15) cm 5.4%	(3.66) cm 3.6%	(3.97) cm 8.4%
15mm	(1.51) cm 8.5%	(2.02) cm 7.9%	(2.51) cm 8.6%	(3.10) cm 7.0%	(3.54) cm 6.6%	(3.94) cm 9.2%
16mm	(1.48) cm 10.6%	(1.98) cm 10.1%	(2.45) cm 10.7%	(3.00) cm 9.8%	(3.48) cm 8.2%	(3.84) cm 11.4%
18mm	(1.43) cm 13.7%	(1.89) cm 14.0%	(2.38) cm 13.3%	(2.88) cm 13.5%	(3.20) cm 15.6%	(3.71) cm 14.5%
20mm	(1.32) cm 20.0%	(1.80) cm 17.9%	(2.30) cm 16.3%	(2.63) cm 20.9%	(2.96) cm 21.9%	(3.27) cm 24.7%
22mm	(1.23) cm 25.4%	(1.67) cm 23.9%	(2.09) cm 24.0%	(2.46) cm 26.0%	(2.85) cm 25.0%	(3.15) cm 27.3%
24mm	(1.05) cm 36.4%	(1.46) cm 33.8%	(1.91) cm 30.3%	(2.07) cm 37.9%	(2.27) cm 40.1%	(2.83) cm 34.9%

CP Stent™ Balloon Sizing Chart

Inner Balloon Pressure (atm)	Stent ID (mm)							
	12mm Diameter RBP = 7.0	14mm Diameter RBP = 6.0	15mm Diameter RBP = 5.0	16mm Diameter RBP = 5.0	18mm Diameter RBP = 4.0	20mm Diameter RBP = 4.0	22mm Diameter RBP = 3.0	24mm Diameter RBP = 3.0
1	2.75	3.22	3.49	3.75	3.94	4.02	4.20	4.28
2	2.85	3.32	3.59	3.85	4.36	4.13	4.33	4.50
3	5.85	6.91	6.89	7.79	8.54	9.20	10.16	10.57
4	6.12	7.00	7.02	7.95	8.71	9.63	10.40	11.08
4.5							10.84	11.94
5	6.20	7.08	7.10	8.04	8.91	10.00		
Outer Balloon Pressure (atm)								
1	10.73	13.08	13.45	14.87	16.85	17.91	20.52	22.79
2	10.86	13.27	14.16	15.10	17.06	18.38	21.46	23.95
3	11.15	13.50	14.55	15.68	17.64	19.42	21.98	24.68
4	11.33	13.68	14.88	15.93	18.06	20.07		
5	11.62	13.87	15.06	16.19				
6	11.80	13.98						
7	12.04							

*This data is based on testing performed using the NuMED BIB® Stent Placement Catheter.

The figures in bold face represent the stent ID @ Rated Burst Pressure.

FOR ALL NUMED CATHETERS AN INFLATION DEVICE WITH PRESSURE GAUGE SHOULD BE USED.

BIB DELIVERY CATHETER BALLOON DIAMETER AND INTRODUCER SIZE	REQUIRED INTRODUCER WITH BARE CP STENT	REQUIRED INTRODUCER WITH COVERED CP STENT
12MM (8F)	10F	12F
14MM (8F)	10F	12F
15MM (9F)	11F	12F
16MM (9F)	11F	12F
18MM (10F)	11F	14F
20MM (10F)	12F	14F
22MM (11F)	12F	14F
24MM (11F)	12F	14F

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G. van Wageningen B.V.